
(Dis)assembly instruction for the GANNET Introducer

REF 10112.09.001

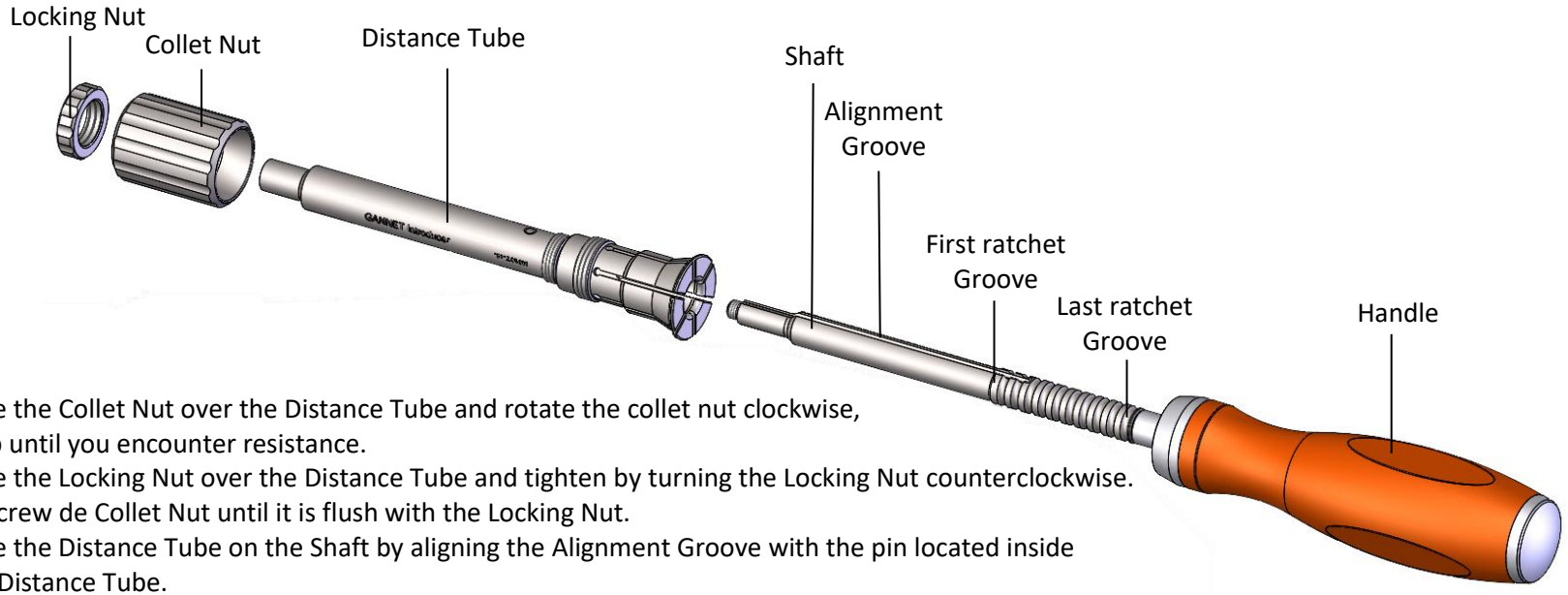

Assembly:

1. Slide the Collet Nut over the Distance Tube and rotate the collet nut clockwise, stop until you encounter resistance.
2. Slide the Locking Nut over the Distance Tube and tighten by turning the Locking Nut counterclockwise.
3. Unscrew de Collet Nut until it is flush with the Locking Nut.
4. Slide the Distance Tube on the Shaft by aligning the Alignment Groove with the pin located inside the Distance Tube.
5. Now tighten the Collet Nut fully over the Collet at the last Ratchet Groove by turning the Collet Nut clockwise.

Disassembly:

1. Untighten the Collet Nut by turning the Collet Nut counterclockwise. The Collet will unlock from the Ratchet Groove.
2. Remove the Distance Tube from the Shaft.
3. Untighten and release the Locking Nut by turning the Locking Nut clockwise.
4. Release the Collet Nut by turning it counterclockwise and removing it from the shaft.

BM-005059-EI-1-Gannet Introducer Dissassembly Instruction • Date of issue: 2025-12-03 • Manufacturer address: Baat Medical Products B.V. • F.Hazemeijerstraat 800 • 7555 RJ Hengelo • The Netherlands • Telephone: +31(0)88 565 66 00 • www.baatmedical.com