
(Dis)assembly instruction for the GANNET Bone Jack

REF 10112.10.001

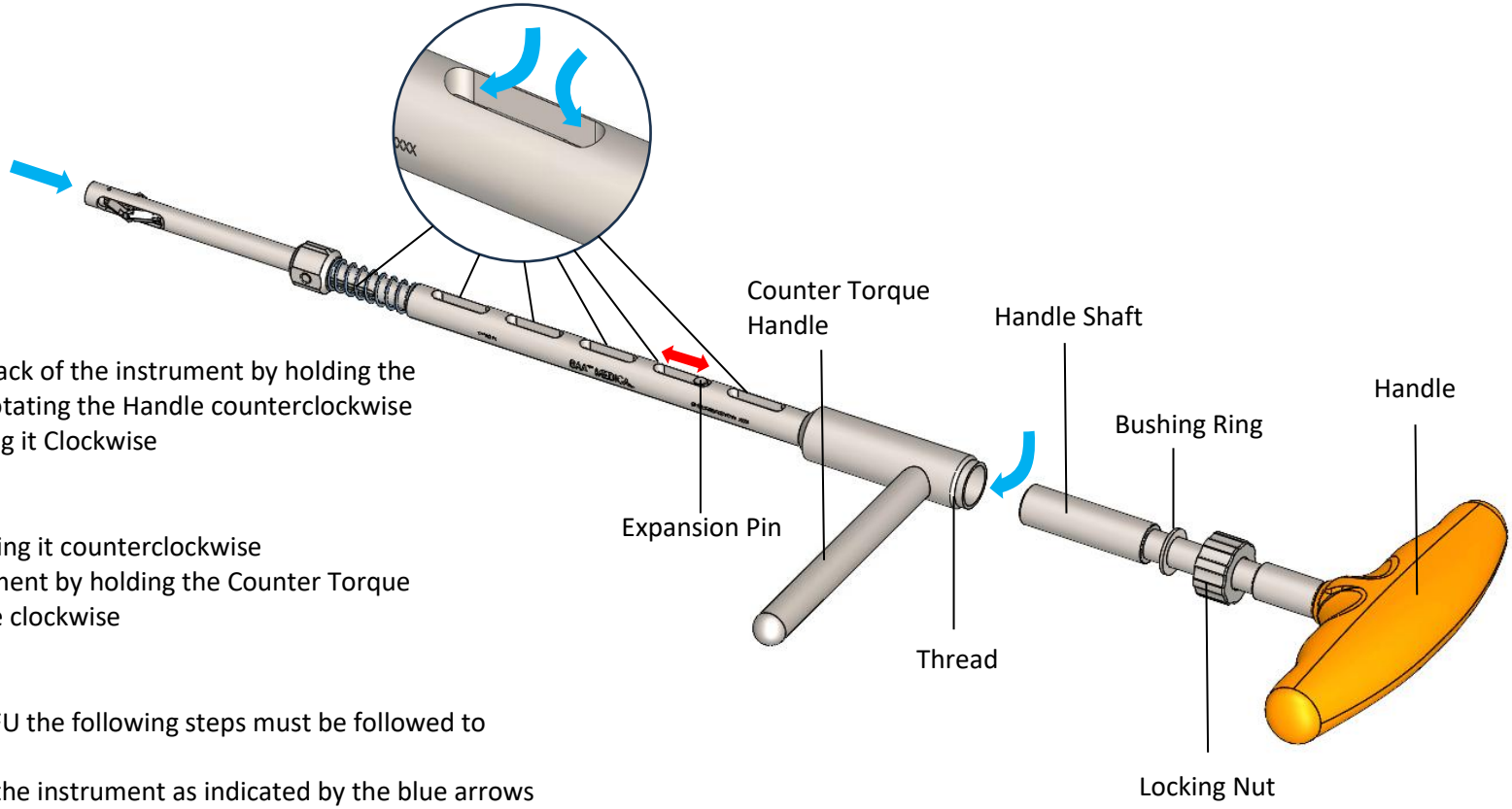

Assembly:

1. Slide the handle shaft in the back of the instrument by holding the Counter Torque Handle and rotating the Handle counterclockwise
2. Lock the Locking Nut by turning it Clockwise

Disassembly:

1. Unlock the Locking Nut by turning it counterclockwise
2. Remove the back of the instrument by holding the Counter Torque Handle and rotating the Handle clockwise

Cleaning

In addition to cleaning instructions in IFU the following steps must be followed to ensure proper cleaning:

1. Use all flush ports to flush out the instrument as indicated by the blue arrows
2. Move Expansion Pin back and forward during flushing indicated by red arrow
3. Rinse thoroughly until all contamination is removed

BM-008336-EI-2-Gannet Bone Jack Disassembly Instructions • Date of issue: 2026-06-01 • Manufacturer address: Baat Medical Products B.V. • F. Hazemeijerstraat 800 • 7555 RJ Hengelo • The Netherlands • Telephone: +31(0)88 565 66 00 • www.baatmedical.com